

Collection of specimens from suspected cases of COVID-19 disease in community settings

Pre-requisite for testing for COVID-19 disease

- All patients being tested for COVID-19 disease **MUST BE NOTIFIED** to the DOH
- Go to Instructions for [Notification of COVID-19](#) or [submit online](#)

Precautions while performing specimen collection

- For most patients with mild illness in the community, collection of upper respiratory specimens is a lower risk procedure and can be performed using **standard, contact and droplet precautions**
- Personal protective equipment (PPE) is to be worn as per the [Department of Health \(DOH\) WA Clinician Alert for General Practitioners](#)

Specimens to collect

Optimal upper respiratory tract specimens are a bilateral **deep nasal** swab (requires flexible shaft Minitip swab) **AND** an **oropharyngeal** swab. If flexible shaft swabs are not available, a rigid shafted PCR swab can be used but this is more uncomfortable for the patient.

Equipment required

- 1 x flexible shaft swab with Universal Transport Media/Virus Transport Media (if available), **OR**
- 1 x dry rigid shaft PCR swab (for mid-turbinate/oropharyngeal swab if flexible shaft swab not available)
- Wooden tongue depressor

Note 1: Serum for COVID-19 serology is not required for patients not admitted to hospital.

How to collect oropharyngeal/throat sample then bilateral deep nasal swab using a single flexible shaft swab

1. Perform hand hygiene before donning PPE.
2. Remove the flexible shaft swab from its packaging.
3. Stand slightly to the side of the patient to avoid exposure to respiratory secretions, should the patient cough or sneeze.
4. **Oropharyngeal/throat sample:** Ask the patient to open their mouth and stick out their tongue. Use a wooden tongue depressor to press the tongue downwards. Firmly swab the posterior pharynx and the tonsillar arches, rotating the swab twice while holding it against the tonsillar fossa at the side of the pharynx. Withdraw the swab, which should be thoroughly wet with throat secretions.
5. **Deep Nasal swab: Using the same swab,** tilt the patient's head back and gently insert the swab into one nostril, then gently insert it along the floor of the nasal cavity parallel to the palate until resistance is encountered (Figure below). Rotate gently for 10-15 seconds, then withdraw, and repeat the process in the other nostril with the same swab, to absorb secretions.
6. Place the swab into the UTM containers, flocked end first, and snap the shaft at the indicated snapping point, using the rim of the UTM tube for leverage if required.
7. Close the UTM container with the flocked end of the swab inside, discarding the remainder of the shaft.
8. Ensure that the UTM container is closed correctly and is not leaking – the swab is usually slightly bent within the tube.

Figure 1: Flocked swab with UTM

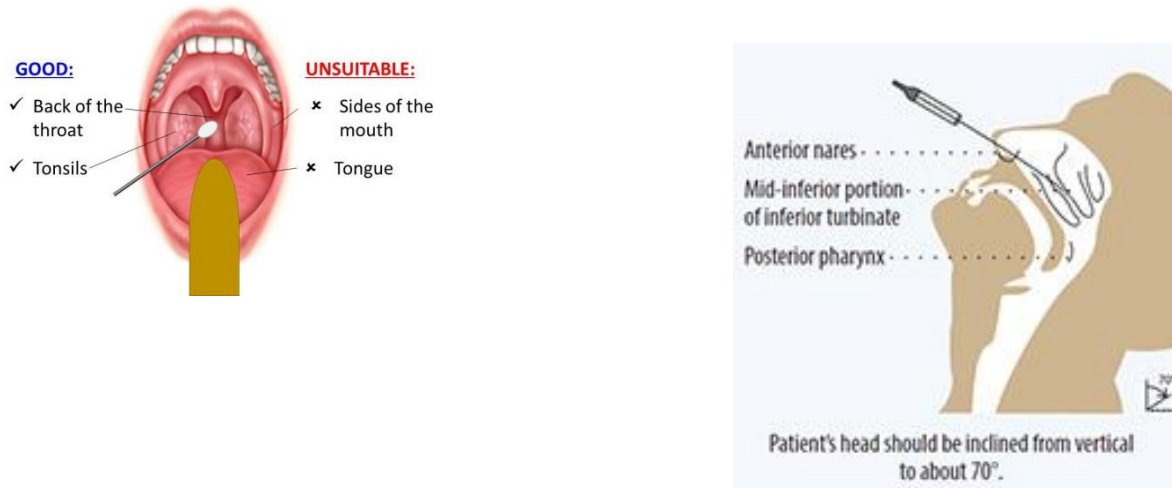


Figure 2: Dry rigid PCR swab



Collection of specimens from suspected cases of COVID-19 disease in community settings

Figure 3: Collection of throat sample swab



Specimen transport

1. Complete a PathWest Request Form with **travel and symptom details, or if the patient is a symptomatic contact of a confirmed case, or if Public Health has requested for the testing to be performed.** Note the samples collected on the form and write 'COVID-19 PCR' on the request form, plus 'respiratory virus PCR' if you require other respiratory virus testing.
2. Label the specimen and place it into the sealed section of a dual compartment biohazard bag. Place the request form in the outer compartment. Do not add any other items to the bag.
3. Send the sample to PathWest as per your usual mechanism.
4. If the testing is considered urgent by yourself or Public Health, call the Clinical Microbiologist for QEII via SCGH switchboard (08 6383 3333).

Figure 4: Collection of deep nasal

Collection of specimens from suspected cases of COVID-19 disease in community settings

Obtaining results

A result is usually available **within 24 hours of receipt of the sample in the PCR laboratory at PathWest on the QEII site**. Call the PathWest Call Centre on 08 6383 4567 for results (do not call the laboratory directly). Please do not ask your patients to call the laboratory for their results as we cannot give results to them.

PLEASE NOTE: COLLECTION OF SAMPLES FROM PATIENTS WITH MODERATE-SEVERE SYMPTOMS

- Please note that it is **NOT** recommended for primary care providers to be collecting lower respiratory tract samples, including sputum samples, from suspected patients with severe symptoms or symptoms suggestive of pneumonia (fever, breathing difficulty, or frequent, severe or productive coughing episodes). It is recommended that these patients should be managed and tested in an Emergency Department.
- See [Coronavirus Disease 2019 \(COVID-19\) CDNA National Guidelines for Public Health Units](#) for instructions for collection in the community setting if referral to hospital for testing of moderate to severely ill patients is not possible. Please note the additional PPE (N95/P2 masks) and other additional precautions are required for this situation.